

# EVACUATION SLYDE™



The Belluscura Evacuation Slyde provides a safe, efficient way to evacuate non-ambulatory occupants from any multistory building using the stairways. It folds flat and stores in stacks of up to five units in a storage sleeve. The sleeve can be wall-mounted for quick access in key locations.



belluscura™

- Folds flat and stacks for easy storage
- Restraints with buckles to secure the occupant for transport
- Corrugated polypropylene, double-wall construction is lightweight and durable
- Strong web handles at corners and sides provide multi-point control for two to seven carriers
- The length of the web handle allows an ergonomically correct upright posture for operators during stair descent

**500LB**  
Capacity

**1YEAR**  
Warranty



# EVACUATION SLYDE™

## Features and Benefits

- Folds flat and stacks 5 deep for easy storage
- Easy set-up and instruction pictograms on each unit
- Four strong web handles for multipoint control by 2 to 7 operators
- Size permits continued stairwell traffic flow during occupant evacuation
- Functions as post-evacuation surge capacity device
- Lightweight for easier transports
- Buckle restraints are quick and easy to use<sup>1</sup>
- Compatible with Belluscura Stairwell Belay System for heavy occupant evacuations.
- Wall-mounted storage options are space-saving and convenient


## Accessories



**Belay Box and Bolt**  
6400-030-000



**Provides stable wall anchor**

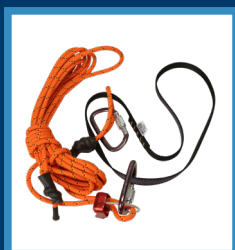

**Stairwell Belay**  
6400-020-000



**5 COUNT**  
6400-010-001

**Slyde Sleeve**

## Specifications

Model Number	6400-000-000
Material	Corrugated polypropylene
Height	8 in. (20.3 cm)
Width <sup>2</sup>	19 in. (48.3 cm)
Length	85 in. (215.9 cm)
Weight	9 lb. (4.1 kg)
Folded Height	2 in. (5.1 cm)
Max. Weight Capacity <sup>3,4</sup>	500 lb. (227 kg)



belluscura™

1. Passenger restraint set includes waist and chest restraints.
2. Dimensions are measured from the outermost edges of the main frame. Specifications are rounded to the nearest whole number. Conversions are calculated before rounding.
3. Operators must consider the weight of the passenger, equipment and accessories when determining the total weight on the product.
4. Maximum weight capacity is total weight distributed in accordance with basic human anatomy. Belluscura reserves the right to change specifications without notice.

## ORDER INFORMATION

email: [customerservice@belluscura.com](mailto:customerservice@belluscura.com)  
tel: +1 214-396-8313  
toll free: +1 844-883-0940

**Belluscura LLC**  
5504 Democracy Drive, Suite 200  
Plano, TX 75024  
United States  
tel: +1 214-396-8313

**Belluscura plc**  
82 St John Street, Farringdon  
London EC1M 4JN  
United Kingdom  
tel: +44 (0) 1865 261506  
toll free: +1 844-883-0940